

A-dec Inspire®
Sterilization Solutions



A-dec Inspire Sterilization Center

Created for dentistry. Made for the dental environment.

Patient protection goes beyond processes and protocols. It's about expediting materials management and streamlining workflow. It considers durable, sustainable materials that withstand wet, harsh dental environments. And it helps ensure compliance by intuitively and safely guiding the dental team through the sterilization process. That's a big ask. But with more than sixty years in the dental industry, thoughtful design based on the way dental teams practice is second nature to A-dec®.



STERILIZATION FRONT AND CENTER

When materials move efficiently, so do you.

The sterilization center is no longer something we hide away in the back room. Along with a pleasing aesthetic, A-dec sterilization centers are organized into functional workspaces, with specified areas for receiving, cleaning, preparation, packaging, sterilization and storage. You can easily group together related steps and simplify the entire sterilization process. The furniture design follows recommended CDC protocol, creating a safe, consistent process and streamlining the flow of materials: Load materials and instruments directly from the sterilization center or sterilized storage, deliver them to the treatment room, and return used materials for processing. Everything is in position and ready when you are, optimizing labor—one of your highest costs.



More than a warranty. It's a commitment.

At A-dec, our focus is on making dentistry better—which includes superior customer service along with superior dental products. Your A-dec operatory equipment is backed by a ten-year warranty, and a lifetime of customer support. We'll be here today, tomorrow, and well into the future. That's our commitment to you.

**For complete details, including exclusions,
visit: [a-dec.com/legal/warranty](https://www.a-dec.com/legal/warranty)**



An efficient, systematic sterilization center helps maintain a safe practice—and saves valuable time.

FOLLOWING PROCESSES IS EASY WHEN THE FLOW IS BU



RECEIVE INSTRUMENTS

An orderly design with a clearly defined entry and exit helps streamline instrument flow and increase efficiency.



PROTECT OPERATORS

Visual cues and separate storage spaces decrease the chance of cross contamination.



MANAGE WASTE

Designated containers safely separate sharps, biowaste and trash, simplifying adherence to safety guidelines.



CLEAN

Designated zones keep dirty instruments and materials from contaminating other areas of the sterilization center.



RINSE

The space to adequately rinse cassettes, baskets or oversized items without splashing is key to this important step.

UILT IN.



LUBRICATE AND FLUSH

A well organized and sequential sterilization center eliminates bottlenecks... and waiting.

PACKAGE

Organized supplies and strategically located supply drawers expedite the packaging process.

STERILIZE

An efficient, systematic sterilization center helps maintain a safe and healthy practice—and saves valuable time.

STORE AND STAGE

Sterilized materials and instrument packages are protected from tears, punctures, water splashes and moisture.

MANAGE MATERIALS

Time-saving organization keeps everything in place and encourages adherence to process and protocol.

WHY CHOOSE A-DEC INSPIRE

Durable materials, made for dental.

Proven quality materials are rigorously tested to withstand the day-to-day stressors of the dental environment. It's hassle-free dependability you can trust.

THE DIFFERENCE IS IN THE DETAILS

Durably and rigorously made. Because what's inside matters.

Everyday cabinets are not designed to withstand the rigors of the dental environment. Sustainably made from top-of-the-line materials, A-dec Inspire is crafted from parts and materials chosen specifically for the dental environment.



High-pressure laminate and thick 2-millimeter, color-matched edge banding on all panels absorbs impact that can happen in the dental environment.



Solid powder-coated steel subbase (rather than stainless cladding) keeps liquids from damaging furniture.



Powder-coated handles won't wear or tarnish over the life of the product.



We bolt our furniture together with metal fasteners, wood screws, bolts, and cams. Using furniture assembling hardware means that if any parts need to be replaced, we can unbolt parts and replace it without having to replace the entire cabinet.



All doors are equipped with best-in-class, soft-close Blum® hinges, adjustable to three closure settings, and cycle tested to ensure 20+ years of use.

PRE-CONFIGURED AND PRACTICAL

Intuitive flow. Practical setup.

Designed with the same durable materials as our A-dec Inspire 500 line, A-dec Inspire 394 offers everything you need to follow CDC guidelines in a simplified, pre-designed unit. You can choose dental cabinetry and countertop colors for a beautiful, blended look.



- Preconfigured setups come in 112" or 126" lengths, with options for ultrasonic sterilization, instrument washer and make-ready setups.
- Ultrasonic cleaning design maximizes production within the flow for small- to medium-sized practices with up to six operatories.
- Upper make-ready storage keeps consumables and sterile instruments ready for the dental operator.



OPTIONS AND ADD-ONS

Add what you need. Get what you want.

Optional add-ons let you target features to increase productivity and optimize workflow—your way. Getting what you want is that easy.

WALL-MOUNTED MATERIALS MANAGEMENT PANEL Create additional customized wall storage with mix and match tip-out bins, hanging bins, and shelving.

COLOR-CODED DRAWER INSERTS Assign a color for each procedure and these tray inserts create a streamlined, stress-free routine.

STERILIZATION TOWER SETUPS Choose the option that fits your requirements: cassette and chamber (slide out tray), dual chamber, or dual chamber and cassette (triple stack).


EXTRA CLEAN/DIRTY INDICATORS Red indicates dirty. Blue indicates clean. Extra visual markers further decrease the chance of cross contamination.

SOAP DISPENSER Add a wall-mounted, motion-activated soap dispenser to eliminate another touchpoint.

OPEN SPACE DESIGN Whether you need a desk for your laptop or cart storage area, open space design offers flexible options.


INSTRUMENT WASHER VERSATILITY Open space design neatly fits a tabletop instrument washer. Or, remove the bottom panel and slide in a full-sized unit.

SINK OPTIONS Choose a single or dual-basin sink, depending on your needs.

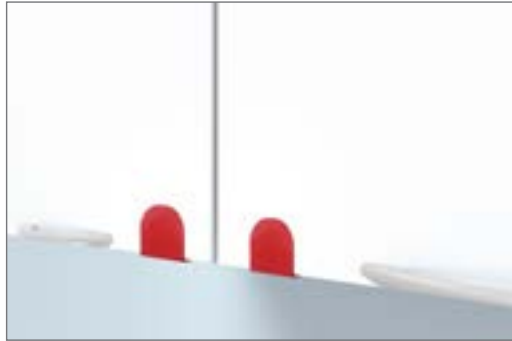

Reliable Creative Solutions

FOR DENTAL PROFESSIONALS WORLDWIDE

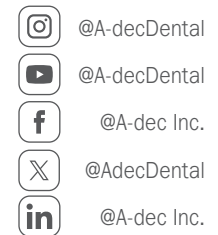

MORE THAN A WARRANTY. IT'S A COMMITMENT.

At A-dec, our focus is on making dentistry better—which includes superior customer service along with superior dental products.

Your A-dec operatory equipment is backed by a ten-year warranty, and a lifetime of customer support. We'll be here today, tomorrow, and well into the future. That's our commitment to you.

For complete details, including exclusions, visit: a-dec.com/legal/warranty

Any colors and products shown may vary and are subject to change without notice. Contact your authorized A-dec dealer for the most current product information. All trademarks in this document are subject to the trademark rights of A-dec, its affiliates, and their collective suppliers in the United States and foreign countries.

©2025 A-dec Inc. All rights reserved. 85.1078.00/PP/1.5M/3-25/Rev A

A-dec Headquarters

2601 Crestview Drive
Newberg, Oregon 97132 USA

P: 1.800.547.1883 within USA/CAN

P: +1.503.538.7478 outside USA/CAN

a-dec.com

Other Locations

A-dec Nashville Experience Center

A-dec Australia

A-dec United Kingdom

